

Computerized pressure Algometer for assessing deep tissue pain perception

Software-based REAL TIME visual and auditory feedback to control and monitor Applied pressure rates

- Real-time visual & auditory feedback
- Available as Stand-alone product or add-on to PATHWAY System
- Comprehensive test statistics
- Body site selection
- USB Patient Response Unit
- Customizable outcome-based color graphic reports
- Calibration verification weight
- Ramp & Hold Protocol



Hardware Advantages


- Ergonomic handle for easy grip
- Rechargeable 50 hour battery
- Support for various tip sizes



Patient Response Unit

Software Capabilities

- Allows managing patients, test paradigms and results
- Pressure units support includes KPa, Kg/cm² and lbs/in²
- Patient response to pressure threshold and tolerance tests can be recorded
- Body site selection
- Test results can be mapped according to a predefined body site
- Test results can be saved, exported to Excel and printed as a customizable report

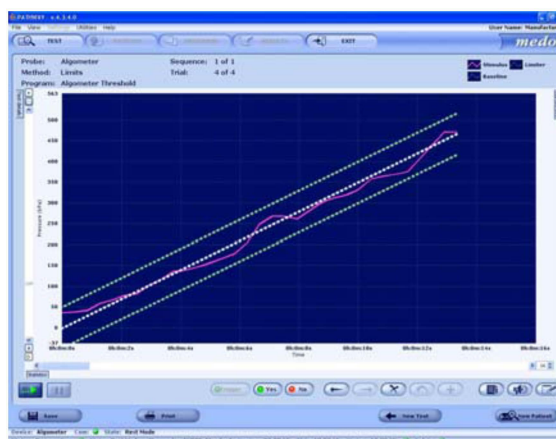

Algometer test screen

Test Management

- Graphic display of test including applied pressure change rate
- Display test statistics according to selected test methods

CE 1282



Refer to operating manual

Medoc Advanced Medical Systems, USA | Compass Medical Technologies, Inc
1502 West Highway 54 - Suite 404, Durham, North Carolina 27707
phone: +972-4-903.8800 | Fax: +972-4-903.8807 | medoc@mindspring.com

Medoc Ltd. Advanced Medical Systems | World Headquarters 1
Ha'Dekel St. Ramat Yishai, Israel 30095
phone: +972-4-903.8800 | Fax: +972-4-903.8808 | info@medoc-web.com

EC REP EU Authorized Representative CEpartner4U BV
Esdoornlaan 13 3951 DB Maan, The Netherlands, office@CEpartner4U.nl
phone: +31-343-442.524 | mobile: +31-6-516-536.26 | Fax: +31-343-442.162

www.medoc-web.com

Medoc
advanced medical systems